



The Next Generation of Crystal Detectors

Ren-Yuan Zhu

California Institute of Technology

May 29, 2015



Why Crystal Calorimeter in HEP?



- Photons and electrons are fundamental particles.
 Precision e/γ measurements enhance physics discovery potential for future HEP experiments.
- Performance of crystal calorimeter in e/γ measurements is well understood:
 - The best possible energy resolution;
 - Good position resolution;
 - Good e/ γ identification and reconstruction efficiency.
- Challenges at future HEP Experiments:
 - Radiation damage at the energy frontier (HL-LHC);
 - Ultra-fast rate at the intensity frontier;
 - Good jet mass resolution at the energy frontier (ILC/CLIC).



Physics with Crystal Calorimeters (II)




Charmonium system observed by CB through Inclusive photons

CB Nal(Tl)

Higgs -> $\gamma\gamma$ by CMS through reconstructing photon pairs



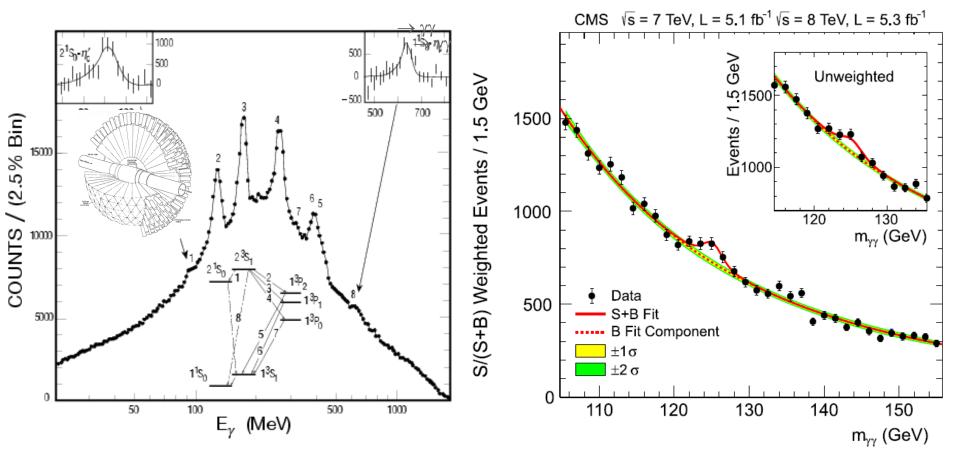




Existing Crystal Calorimeters



Date	75-85	80-00	80-00	80-00	90-10	94-10	94-10	95-20
Experiment	C. Ball	L3	CLEO II	C. Barrel	KTeV	BaBar	BELLE	CMS
Accelerator	SPEAR	LEP	CESR	LEAR	FNAL	SLAC	KEK	CERN
Crystal Type	Nal(TI)	BGO	CsI(TI)	CsI(TI)	Csl	CsI(Tl)	CsI(Tl)	PbWO ₄
B-Field (T)	-	0.5	1.5	1.5	-	1.5	1.0	4.0
r _{inner} (m)	0.254	0.55	1.0	0.27	-	1.0	1.25	1.29
Number of Crystals	672	11,400	7,800	1,400	3,300	6,580	8,800	76,000
Crystal Depth (X ₀)	16	22	16	16	27	16 to 17.5	16.2	25
Crystal Volume (m ³)	1	1.5	7	1	2	5.9	9.5	11
Light Output (p.e./MeV)	350	1,400	5,000	2,000	40	5,000	5,000	2
Photosensor	PMT	Si PD	Si PD	WS ^a +Si PD	PMT	Si PD	Si PD	APD^a
Gain of Photosensor	Large	1	1	1	4,000	1	1	50
σ_N /Channel (MeV)	0.05	0.8	0.5	0.2	small	0.15	0.2	40
Dynamic Range	104	10 ⁵	104	104	104	10 ⁴	10 ⁴	10 ⁵

Future Crystal Calorimeters in HEP: LYSO for COMET (Mu2e, Super B) and CMS at HL-LHC BaF₂ and PbF₂ for Mu2e and g-2 respectively at Fermilab PbF₂, PbFCl and BSO for Homogeneous HCAL for ILC

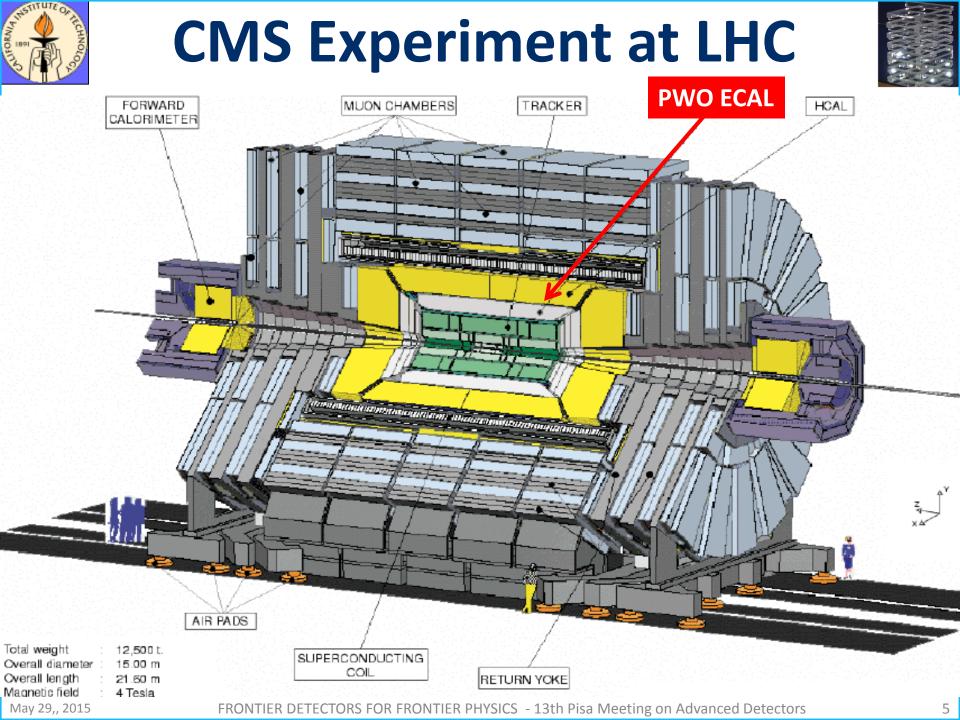




relative response

CMS PWO Monitoring Response



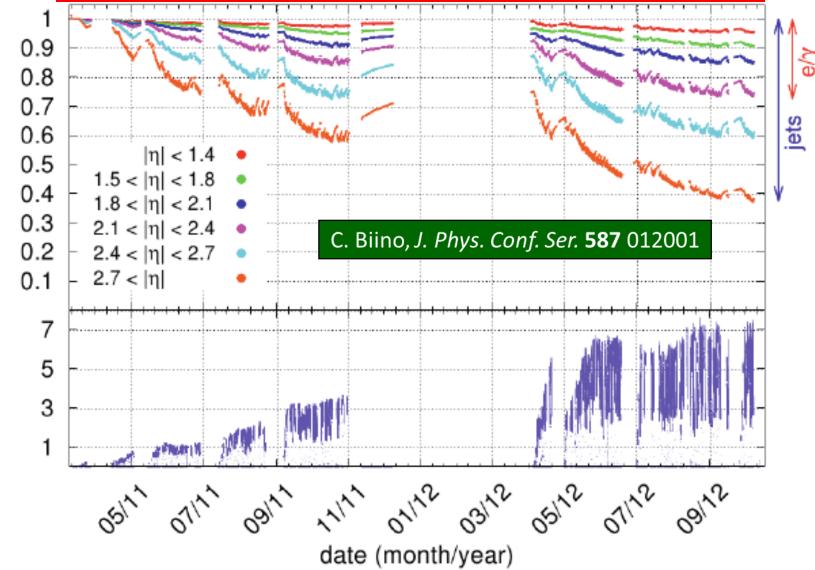

May 29,, 2015

L (10³³ cm⁻² s⁻¹)



Dose Rate Dependent Damage in LO



IEEE Trans. Nucl. Sci., Vol. 44 (1997) 458-476


Light output reaches an equilibrium during irradiations under a defined dose rate, showing dose rate dependent radiation damage

$$dD = \sum_{i=1}^{n} \{-a_i D_i dt + (D_i^{all} - D_i) b_i R dt\}$$

$$D = \sum_{i=1}^{n} \{ \frac{b_i R D_i^{all}}{a_i + b_i R} \left[1 - e^{-(a_i + b_i R)t} \right] + D_i^0 e^{-(a_i + b_i R)t} \}$$

- D_i : color center density in units of m⁻¹;
- D_i^0 : initial color center density;
- D_i^{all} is the total density of trap related to the color center in the crystal;
- a_i : recovery costant in units of hr⁻¹;
- b_i : damage contant in units of kRad⁻¹;
- R: the radiation dose rate in units of kRad/hr.

$$D_{eq} = \sum_{i=1}^{n} \frac{b_i R D_i^{all}}{a_i + b_i R}$$

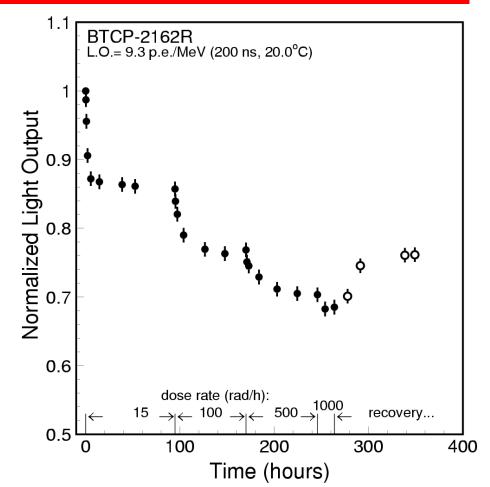


May 29,, 2015



Oxygen Vacancies Identified by TEM/EDS

5 to10 nm black spots identified by TOPCON-002B scope, 200 kV, 10 uA Localized stoichiometry analysis by JEOL JEM-2010 scope and Link ISIS EDS

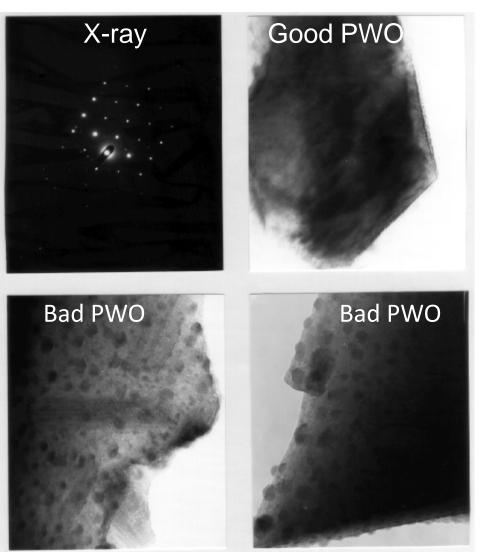

NIM A413 (1998) 297

Atomic Fraction (%) in PbWO₄

As Grown Sample

Element	Black Spot	Peripheral	$Matrix_1$	Matrix ₂
0	1.5	15.8	60.8	63.2
W	50.8	44.3	19.6	18.4
Pb	47.7	39.9	19.6	18.4

The Same Sample after Oxygen Compensation


Element	$Point_1$	$Point_2$	Point ₃	Point ₄
0	59.0	66.4	57.4	66.7
W	21.0	16.5	21.3	16.8
Pb	20.0	17.1	21.3	16.5

May 29,, 2015



Prediction of PWO Radiation Damage



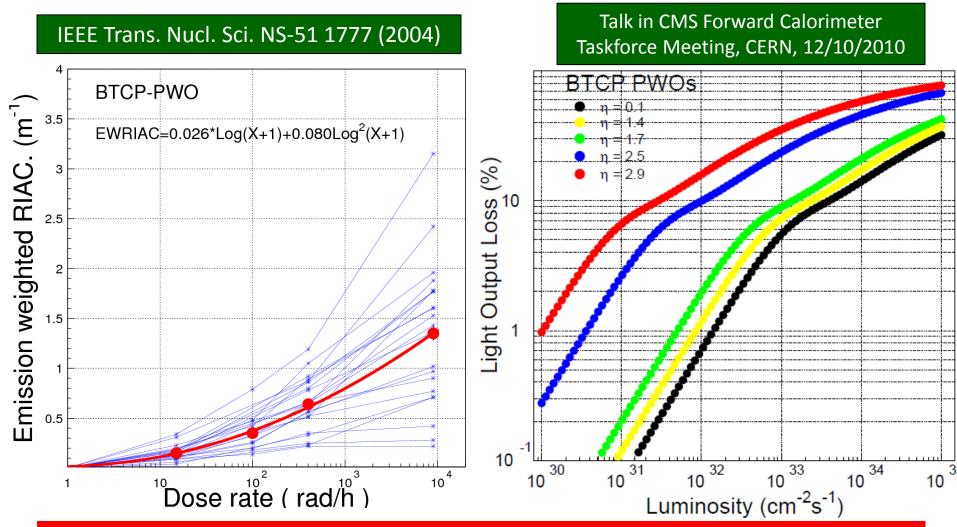


Predicted EM dose induced damage agrees well with the LHC data In addition, there is cumulative hadron induced damage in PWO

May 29,, 2015



Proton Induced Damage in PWO



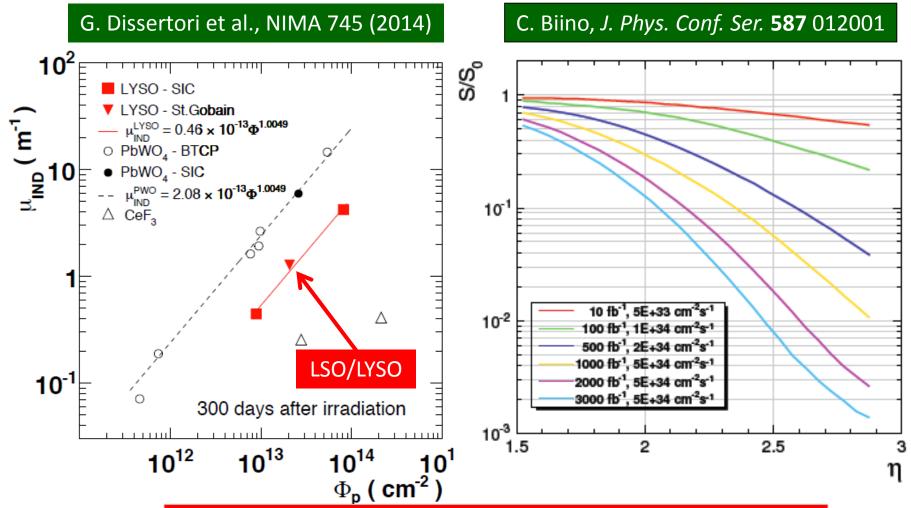


Expected LO loss would be very large for PWO Proton induced absorption in LYSO is 1/5 of PWO



Hadron Fluence @ 3,000 fb⁻¹



Expected neutral and charged hadron fluence is 5×10^{15} /cm² and 3×10^{14} /cm² respectively for the proposed Shashlik endcap at $|\eta| = 3$

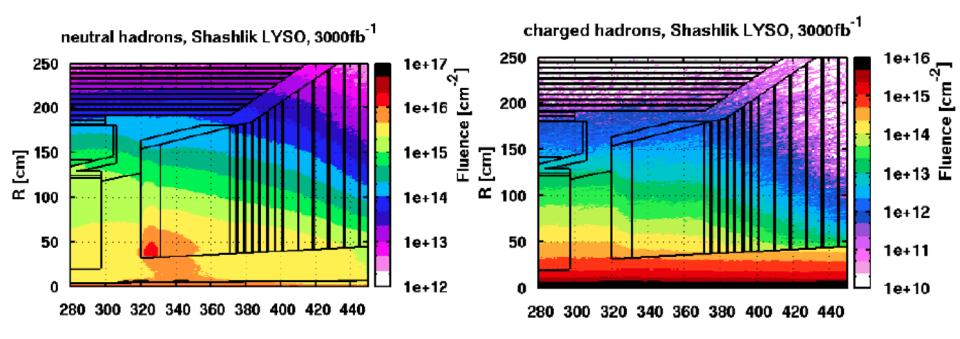

No experimental data show that hadrons (charged or neutral) of > 20 MeV would damage scintillators equally, so they are treated separately



Expected Fluence from FLUKA



CMS Radiation	LHC (10 ³⁴ cm ⁻¹	² s ⁻¹ , 500 fb ⁻¹)	HL-LHC (5×10 ³⁴ cm ⁻² s ⁻¹ , 3000 fb ⁻¹)		
Civio Naciación	Barrel (max)	Endcap (max)	Barrel (max)	Endcap (max)	
Absorbed dose (rad)	3.50E+05	2.10E+07	2.10E+06	1.26E+08	
Dose rate (rad/h)	25	1512	126	7560	
Fast neutrons fluence (E>100KeV, cm ⁻²)	3.00E+13	8.00E+14	1.80E+14	4.80E+15	
Fast neutrons flux (E>100KeV, cm ⁻² s ⁻¹)	6.00E+05	1.60E+07	3.00E+06	8.00E+07	
Charged hadrons fluence (cm ⁻²)	4.00E+11	5.00E+13	2.40E+12	3.00E+14	
Charged hadrons flux (cm ⁻² s ⁻¹)	8.00E+03	1.00E+06	4.00E+04	5.00E+06	


g-rays: Up to 130 Mrad at 7.6 krad/h, See [1] and [2] Fast Neutrons: Up to 5 x 10^{15} n/cm² at 8 x 10^7 n/cm²/s Charged hadrons: Up to 3 x 10^{14} p/cm² at 5 x 10^6 p/cm²/s



Energy Spectra Expected at HL-LHC



The peak energy of charged hadrons at CMS endcap is hundreds MeV Proton energy used for irradiations: 67 MeV/800 MeV/24GeV @ Davis/LANL/CERN Neutron energy used in irradiation: 2.5 MeV from a pair of Cf-252 source

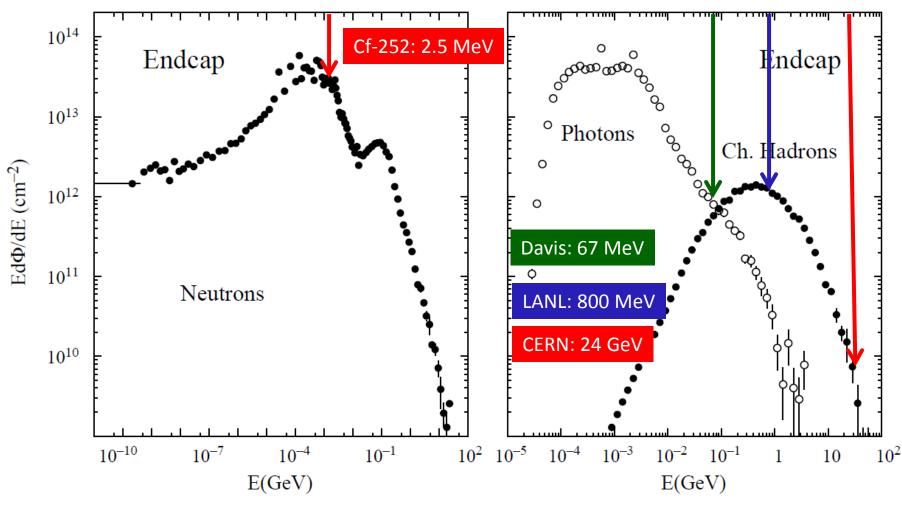

FRONTIER DETECTORS FOR FRONTIER PHYSICS - 13th Pisa Meeting on Advanced Detectors



No Neutron Damage in PWO



5.2 Radiation damage effects under neutron irradiation


In view of the intense neutron flux expected in CMS (see section 2) the effects on lead tungstate of neutron exposure were studied in nuclear reactors [47, 48]. The neutron fluxes and energies in these exposures were comparable to those expected in CMS. However, in reactors there is a strong associated gamma dose. The effect arising from neutrons was estimated by comparing the reactor results with results obtained from pure gamma irradiations. This indicated that there was no specific effect due to neutrons on the optical and scintillating properties of lead tungstate, at least up to fluences of 10^{14} cm⁻². This was confirmed by later independent studies [49]. It is also to be mentioned that recent tests performed at a very high fluence, of the order of 10^{19} to 10^{20} n·cm⁻² and 330 MGy (i.e. well above the level that will be ever achieved in any physics experiment) revealed the robustness of lead tungstate crystals which were not destroyed nor locally vitrified, and remained scintillating after such heavy irradiation [50].

The CMS Electromagnetic Calorimeter Group, *Radiation hardness qualification of PbWO*₄ scintillation crystals for the CMS Electromagnetic Calorimeter, 2010 JINST 5 P03010



A Saclay Paper on Neutron Damage to 10¹⁹ n/cm²



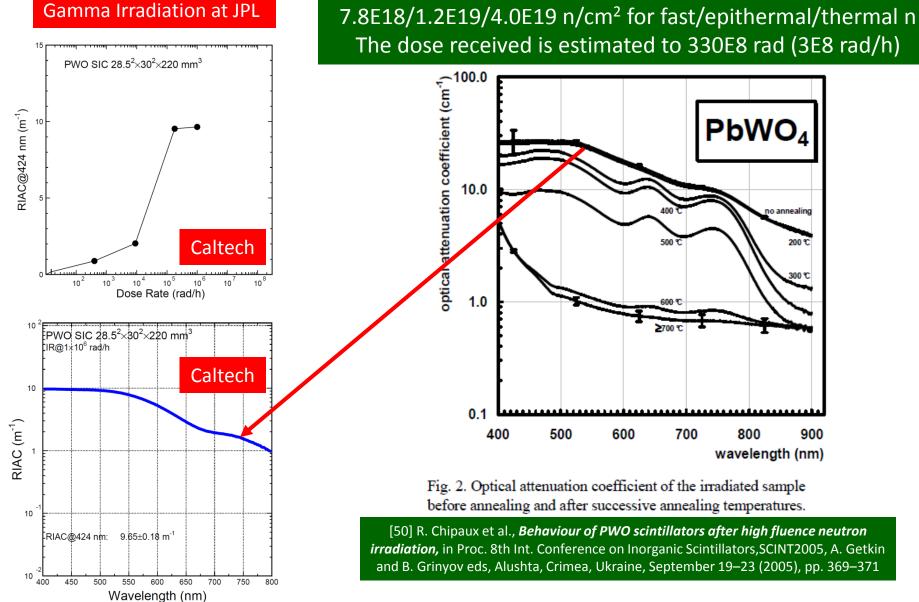

May 29,, 2015



Bright, Fast Scintillator: LSO/LYSO



	Crystal	Nal(TI)	CsI(TI)	Csl	BaF ₂	BGO	LYSO(Ce)	PWO	PbF ₂
	Density (g/cm ³)	3.67	4.51	4.51	4.89	7.13	7.40	8.3	7.77
	Melting Point (°C)	651	621	621	1280	1050	2050	1123	824
	Radiation Length (cm)	2.59	1.86	1.86	2.03	1.12	1.14	0.89	0.93
	Molière Radius (cm)	4.13	3.57	3.57	3.10	2.23	2.07	2.00	2.21
	Interaction Length (cm)	42.9	39.3	39.3	30.7	22.8	20.9	20.7	21.0
	Refractive Index ^a	1.85	1.79	1.95	1.50	2.15	1.82	2.20	1.82
	Hygroscopicity	Yes	Slight	Slight	No	No	No	No	No
	Luminescence ^b (nm) (at peak)	410	550	310	300 220	480	402	425 420	?
	Decay Time ^b (ns)	245	1220	26	650 0.9	300	40	30 10	?
	Light Yield ^{b,c} (%)	100	165	3.7	36 4.1	21	85	0.3 0.1	?
	d(LY)/dT ^ь (%/ ºC)	-0.2	0.4	-1.4	-1.9 0.1	-0.9	-0.2	-2.5	?
	Experiment	Crystal Ball	BaBar BELLE BES III	KTeV	(GEM) TAPS Mu2e	L3 BELLE EIC?	Comet {Mu2e,SuperB) CMS?	CMS ALICE PANDA	A4 g-2 HHCAL?
	a. at peak of emis	sion; b. up/	low row: sl	ow/fast con	nponent; d		dout device tak		
Ma	ay 29,, 2015 FR	ONTIER DETEC	TORS FOR FR	ONTIER PHYSI	CS - 13th Pi	sa Meeting o	n Advanced Detect	tors	16



Bright, Fast & Rad Hard LYSO





623 SIC LYSO Plates of 14 x 14 x 1.5 mm with Five Holes

- LYSO is a bright (200 times of PWO), fast (40 ns) and radiation hard crystal scintillator.
- No variation in emission spectrum was observed after y-rays irradiation, indicating scintillation mechanism is not damaged.
- y-ray induced radiation damage in LYSO does not recover at room temperature, indicating a stable calorimeter *in situ*.
- γ-ray induced absorption coefficient measured for 20 cm long LYSO crystals after 200 Mrad irradiation is about 2 m⁻¹.
- γ-ray induced light output loss measured for 14 x 14 x 1.5 mm plates after 100 and 200 Mrad irradiation is about 6 and 8% respectively.
- The material is widely used in the medical industry with existing mass production capability.




y-Ray Induced Damage in LYSO



No damage in scintillation

No recovery in room temperature

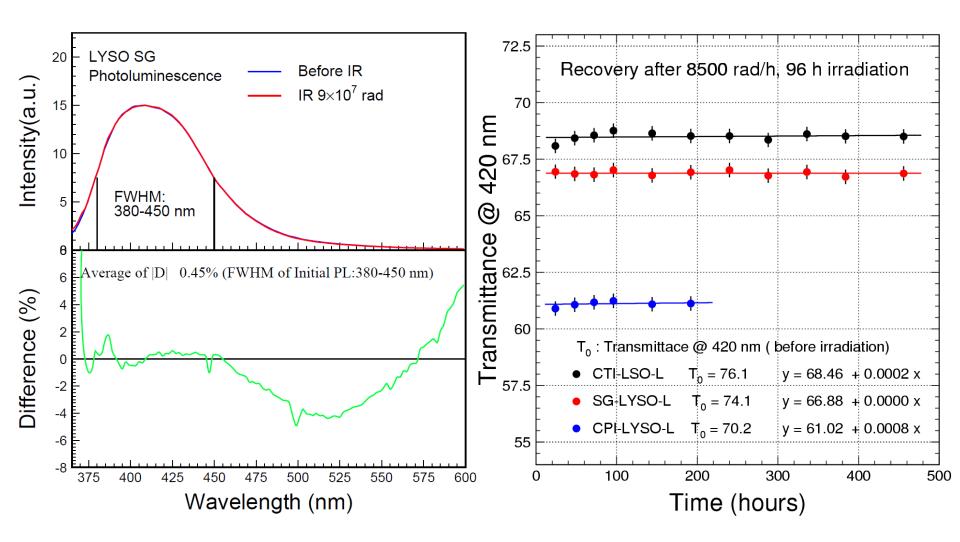




Radiation Hardness of 20 cm LYSO




About 20% loss and <5% divergence for six vendors after 10 Mrad

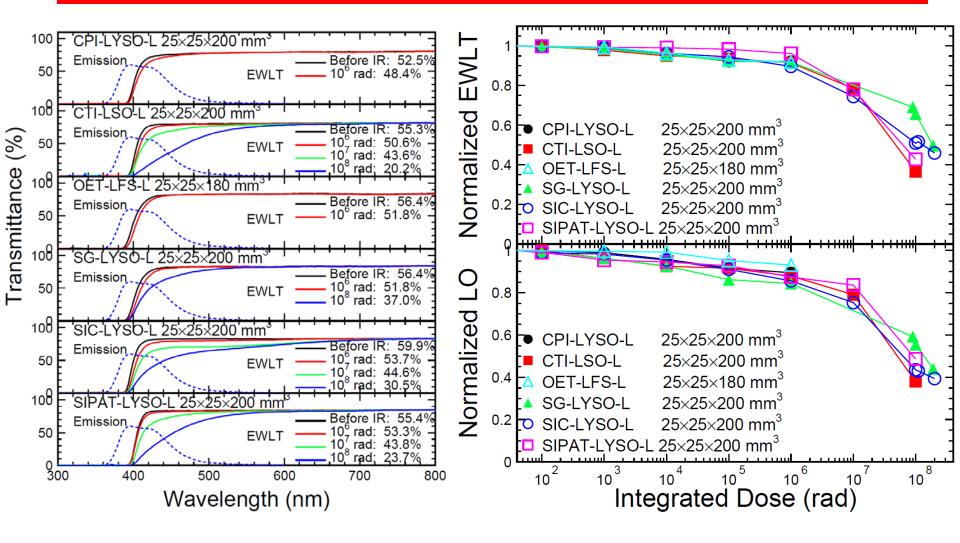



An Option for CMS FCAL Upgrade



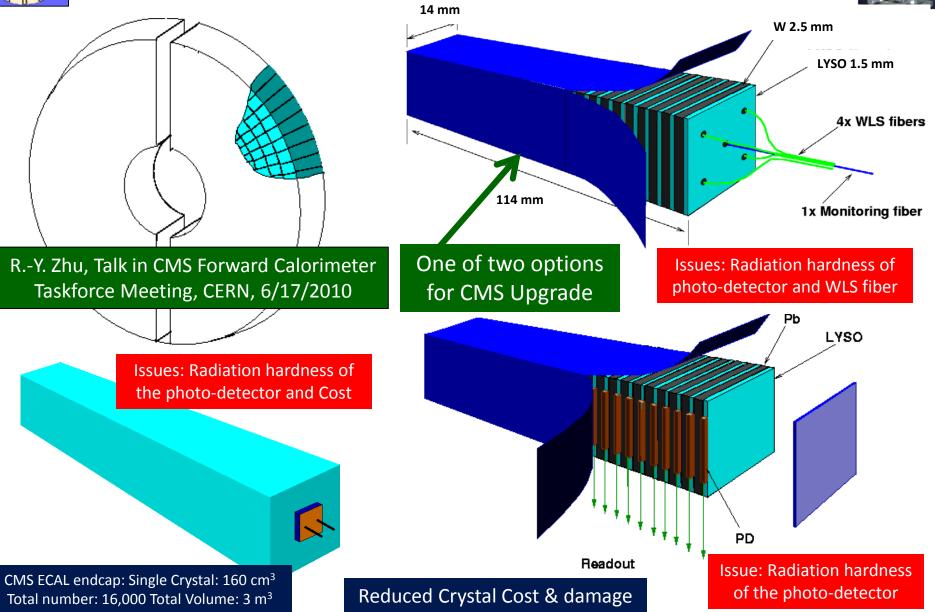



May 29,, 2015



A Shashlik Cell Irradiated at JPL



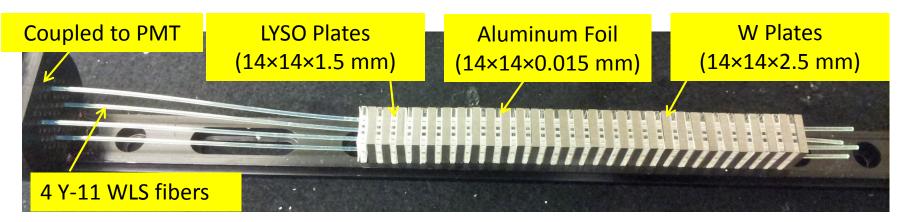

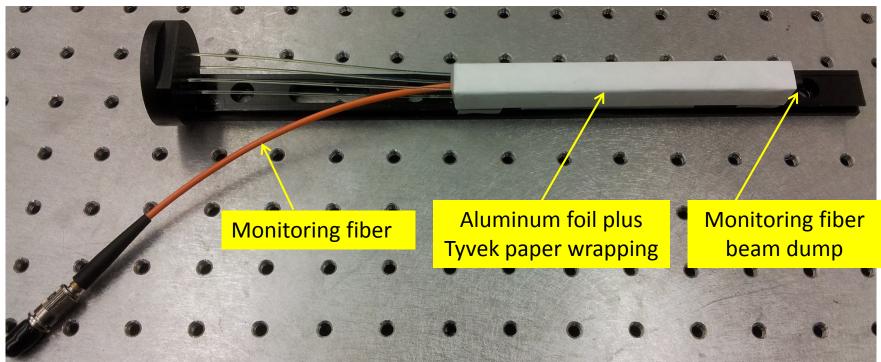


May 29,, 2015





EWRIAC = 1.5, 3 and 4 m⁻¹ after 10, 100 and 180 Mrad LO loss after 100 Mrad is 4 and 6% respectively for direct and WLS readout

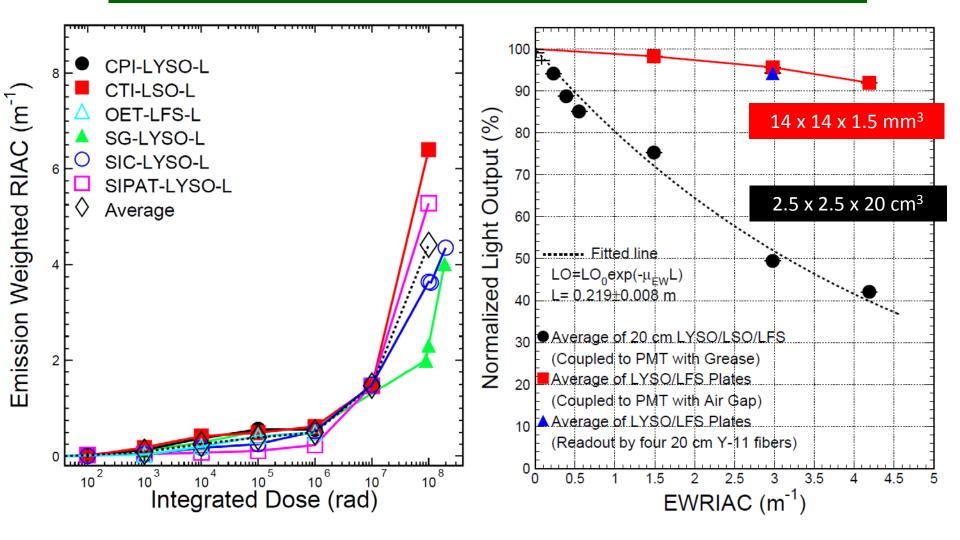




Summary of Proton Induced Damage




A 20 cm and four 14×14×1.5 mm³ LYSO were irradiated by 800 MeV and 24 GeV protons at LANL and CERN respectively. The result shows that the expected RIAC at the HL-LHC is about a few m⁻¹, indicating loss of 4 and 6% respectively for direct and WLS readout.

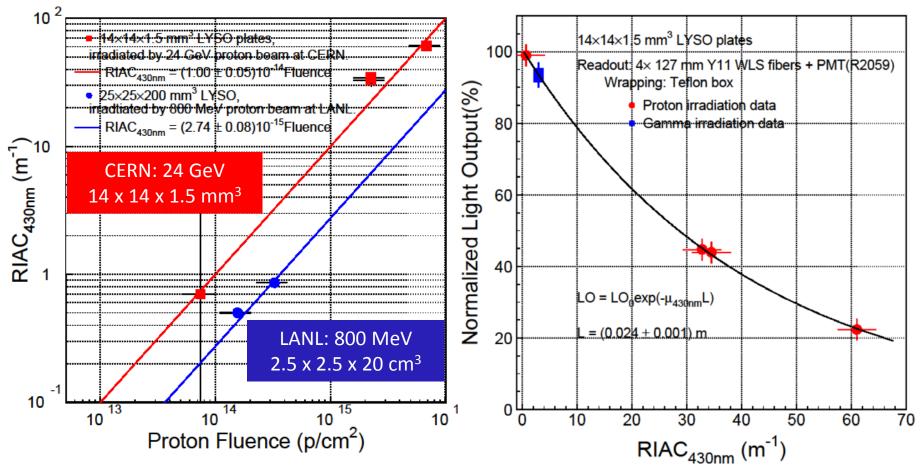



LSO/LYSO Crystal Cost



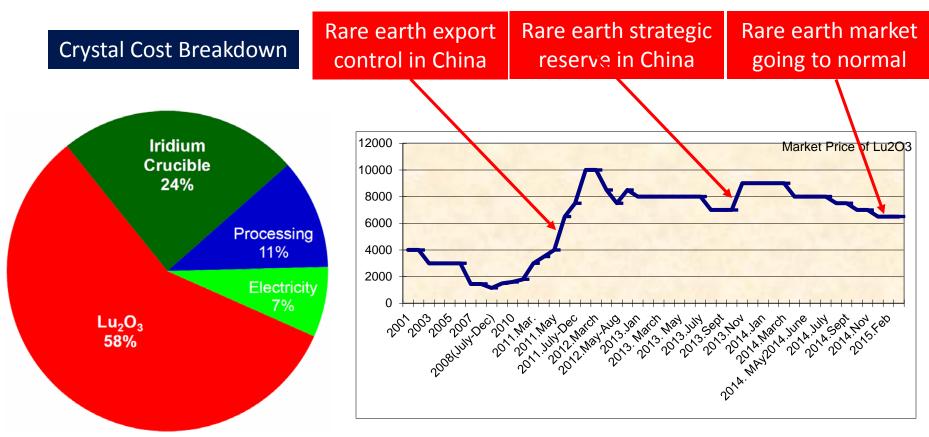

Assuming Lu₂O₃ at \$400/kg and 33% yield the cost is about \$18/cc. Quotations received at \$22-25/cc.

Current Lu_2O_3 price indicates that LYSO price is going down from \$42/cc last year



Alternative Fast Crystals



Talk in CMS Forward Calorimetry Task Force Meeting, CERN, June 27, 2012

	LSO/LYSO	GSO	YSO ¹	Csl	BaF ₂	CeF ₃	CeBr ₃ 2	LaCl ₃	LaBr ₃	Plastic scintillator (BC 404) [€]
Density (g/cm ³)	7.40	6.71	4.44	4.51	4.89	6.16	5.23	3.86	5.29	1.03
Melting point (°C)	2050	1950	1980	621	1280	1460	722	858	783	70#
Radiation Length (cm)	1.14	1.38	3.11	1.86	2.03	1.70	1.96	2.81	1.88	42.54
Molière Radius (cm)	2.07	2.23	2.93	3.57	3.10	2.41	2.97	3.71	2.85	9.59
Interaction Length (cm)	20.9	22.2	27.9	39.3	30.7	23.2	31.5	37.6	30.4	78.8
Z value	64.8	57.9	33.3	54.0	51.6	50.8	45.6	47.3	45.6	-
dE/dX (MeV/cm)	9.55	8.88	6.56	5.56	6.52	8.42	6.65	5.27	6.90	2.02
Emission Peak ^a (nm)	420	430	420	420 310	300 220	340 300	371	335	356	408
Refractive Index ^b	1.82	1.85	1.80	1.95	1.50	1.62	1.9	1.9	1.9	1.58
Relative Light Yield ^{a,c}	100	45	76	4.2 1.3	42 4.8	8.6	141	15 49	153	35
Decay Time ^a (ns)	40	73	60	30 6	650 0.9	30	17	570 24	20	1.8
d(LY)/dT ^d (%/°C)	-0.2	-0.4	-0.3	-1.4	-1.9 0.1	~0	-0.1	0.1	0.2	~0

a.

- At the wavelength of the emission maximum. b.
- Top line: slow component, bottom line: fast component.¹. N. Tsuchida et al Nucl. Instrum. Methods Phys. Res. A, 385 (1997) 290-298 http://www.hitachi-chem.co.jp/english/products/cc/017.html

2. W. Drozdowski et al. IEEE TRANS. NUCL. SCI, VOL.55, NO.3 (2008) 1391-1396 Chenliang Li et al, Solid State Commun, Volume 144, Issues 5-6 (2007),220-224 http://scintillator.lbl.gov/

- 3. http://www.detectors.saint-gobain.com/Plastic-Scintillator.aspx http://pdg.lbl.gov/2008/AtomicNuclearProperties/HTML PAGES/216.html
- Relative light yield normalized to the light yield of LSO c.

d. At room temperature (20°C)

Softening point

Rising Time for 1.5 X₀ **Samples**



Talk in the time resolution workshop at U. Chicago, 4/28/2011: Agilent MSO9254A (2.5 GHz) DSO with 0.14 ns rise time Hamamatsu R2059 PMT (2500 V) with rise time 1.3 ns

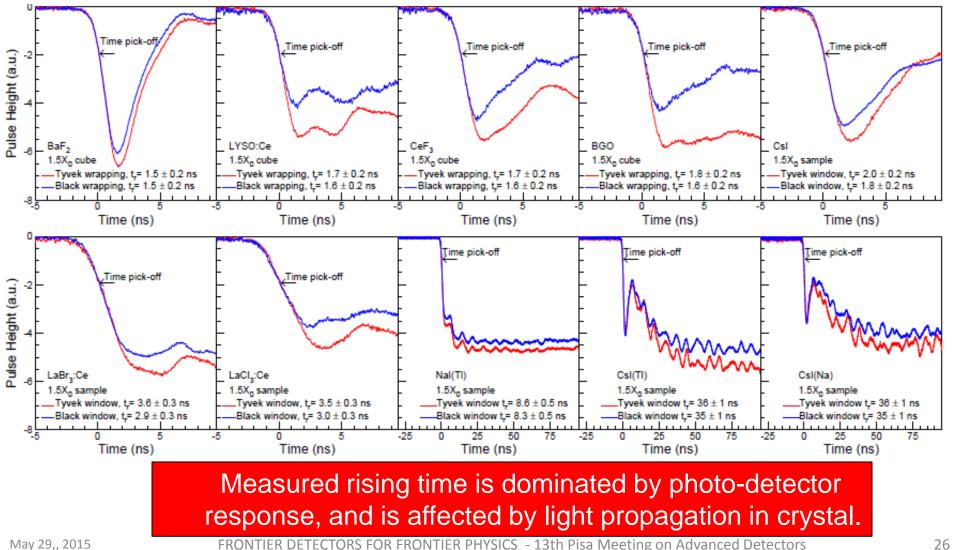




Figure of Merit for Timing



FoM is calculated as the LY in 1^{st} ns obtained by using light output and decay time data measured for 1.5 X₀ crystal samples.


Crystal Scintillators	Relative LY (%)	A ₁ (%)	τ ₁ (ns)	A ₂ (%)	τ ₂ (ns)	Total LO (p.e./MeV, XP2254B)	LO in 1ns (p.e./MeV, XP2254B)	LO in 0.1ns (p.e./MeV, XP2254B)	LY in 0.1ns (photons/MeV)
BaF ₂	40.1	91	650	9	0.9	1149	71.0	11.0	136.6
LSO:Ca,Ce	94	100	30			2400	78.7	8.0	110.9
LSO/LYSO:Ce	85	100	40			2180	53.8	5.4	75.3
CeF ₃	7.3	100	30			208	6.8	0.7	8.6
BGO	21	100	300			350	1.2	0.1	2.5
PWO	0.377	80	30	20	10	9.2	0.42	0.04	0.4
LaBr ₃ :Ce	130	100	20			3810	185.8	19.0	229.9
LaCl ₃ :Ce	55	24	570	76	24	1570	49.36	5.03	62.5
Nal:Tl	100	100	245			2604	10.6	1.1	14.5
Csl	4.7	77	30	23	6	131	7.9	0.8	10.6
CsI:TI	165	100	1220			2093	1.7	0.2	4.8
Csl:Na	88	100	690			2274	3.3	0.3	4.5

The best crystal scintillator for ultra-fast timing is BaF₂ and LSO(Ce/Ca) and LYSO(Ce). LaBr₃ is a material with high potential.



Mu2e BaF₂ Calorimeter



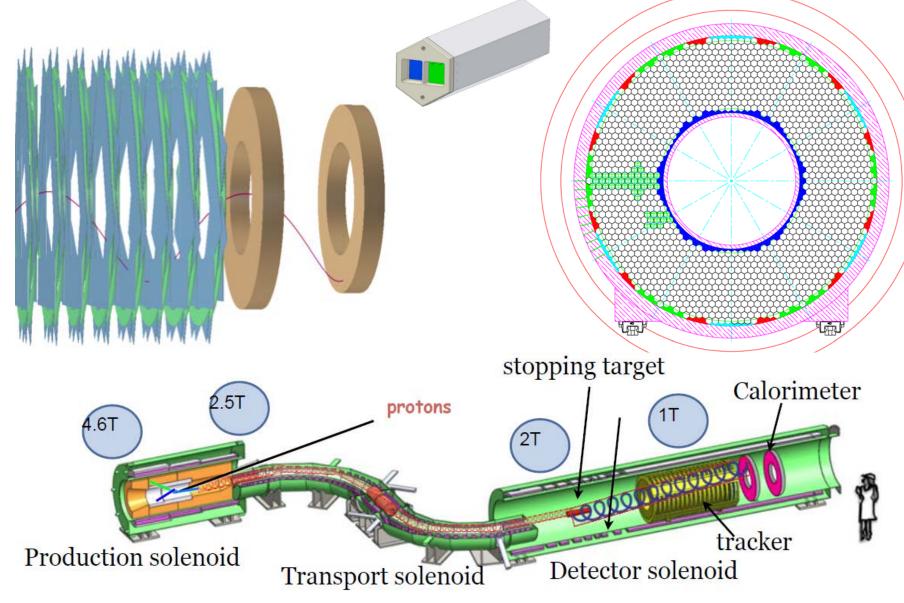




BaF₂ for Very Fast Calorimeter

The Light output of the fast component of BaF₂ crystals at 220 nm with sub-ns decay time is similar to pure CsI.

Spectroscopic selection of fast component may be achieved with solar blind photocathode and/or selective doping.

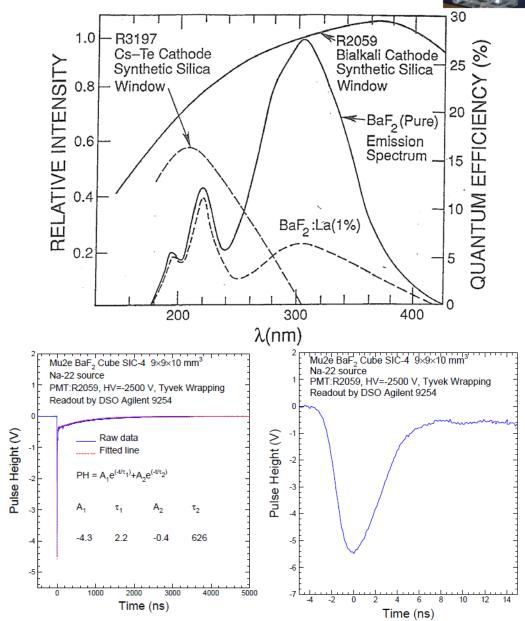




Slow Suppression by Doping and Readout



Y or La doping is effective in improving the F/S ratio for Ba_{0.9}R_{0.1}F₂ powders

B.P. SOBOLEV et al., "SUPPRESSION OF BaF2 SLOW COMPONENT OF X-RAY LUMINESCENCE IN NON-STOICHIOMETRIC Ba0.9R0.1F2 CRYSTALS (R=RARE EARTH ELEMENT)," Proceedings of The Material Research Society: Scintillator and Phosphor Materials, pp. 277-283, 1994.

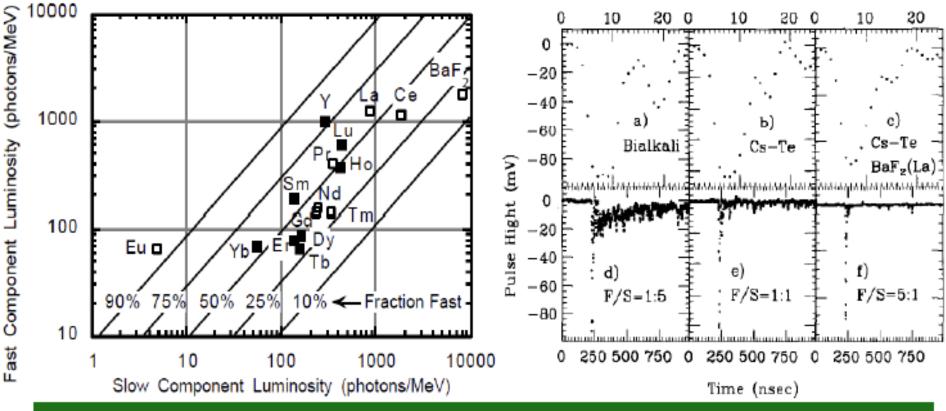

Solar blind cathode is also effective. R&D on doping will be carried out in 2015.


Z. Y. Wei, R. Y. Zhu, H. Newman, and Z. W. Yin, "Light Yield and Surface-Treatment of Barium Fluoride-Crystals," Nucl Instrum Meth B, vol. 61, pp. 61-66, Jul 1991.



Delta-doping for CCD detectors



D. Hitlin *et al*, in this meeting

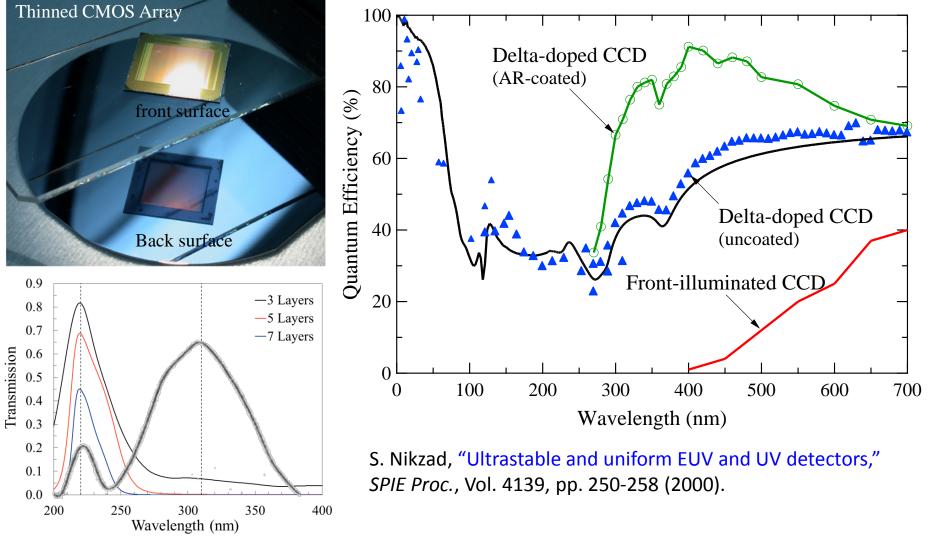




Damage in Long BaF₂ Crystals



Radiation damage in BaF2 crystals saturates at a few tens of krad

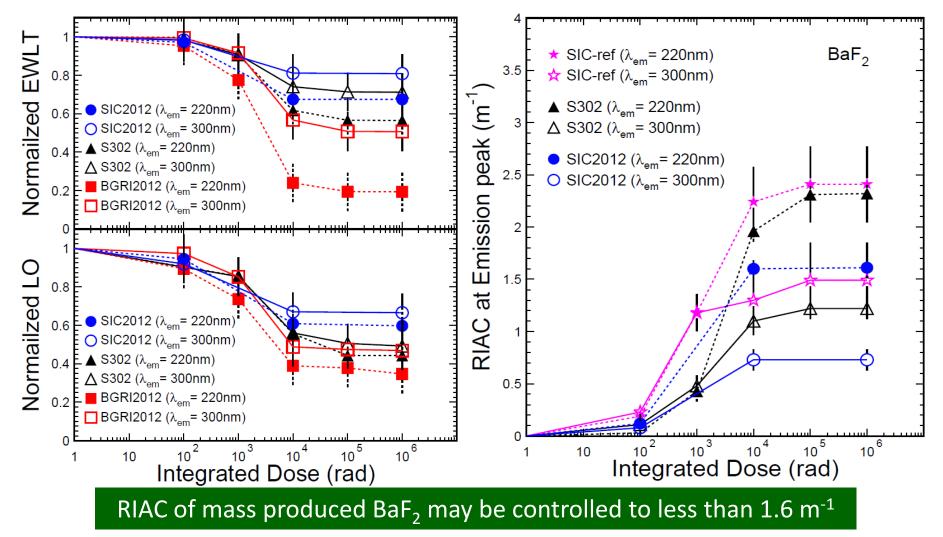


May 29,, 2015



Damage in Long Pure Csl Crystals



Consistent damage between 30/20 cm long pure CsI from SIC/Kharkov

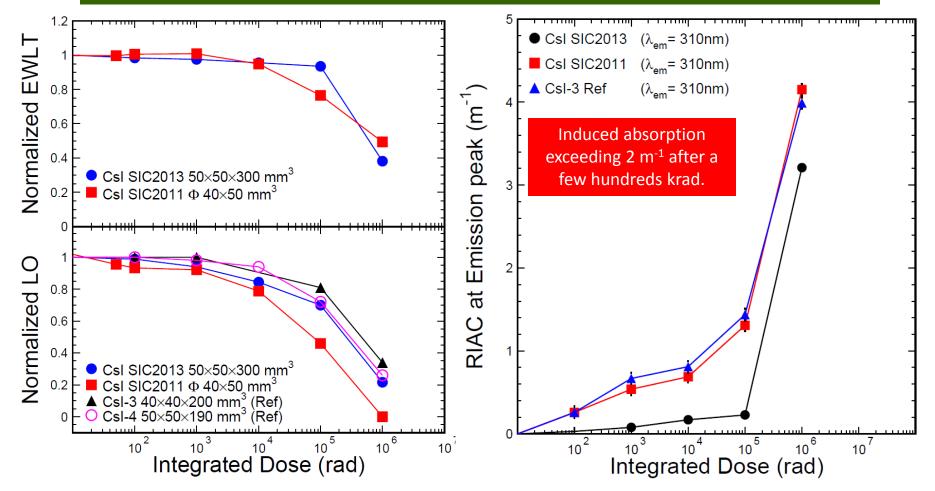


Data of Kharkov crystals: Nucl. Ins. Meth. A 326 (1993) 508-512



Homogeneous Hadronic Calorimeter



A Fermilab team (A. Para et al.) proposed a total absorption homogeneous hadronic calorimeter (HHCAL) detector concept to achieve good jet mass resolution by measuring both Cherenkov and Scintillation light.

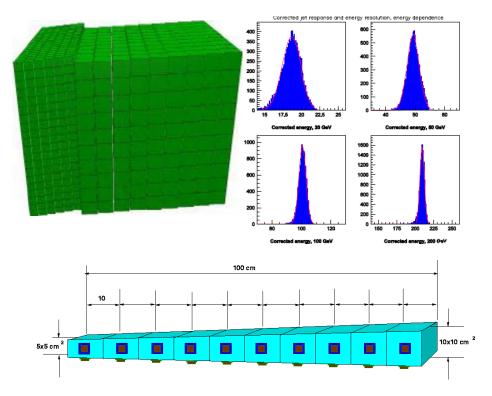

ILCWS-08, Chicago: a HHCAL cell with pointing geometry

Requirements for the Materials:

- Cost-effective material: for 70~100 m³
- Short nuclear interaction length: ~ 20 cm.
- Good UV transmittance: UV cut-off < 350 nm, for readout of Cherenkov light.
- Some scintillation light, not necessary bright and fast.
- Discrimination between Cherenkov and scintillation lights, in spectral or temporal domain.



Candidate Crystals for HHCAL




Cost-effective, UV transparent crystals with both scintillation and Cherenkov light


Parameters	Bi ₄ Ge ₃ O ₁₂ (BGO)	PbWO ₄ (PWO)	PbF ₂	PbClF	Bi ₄ Si ₃ O ₁₂ (BSO)			
ρ (g/cm³)	7.13	8.29	7.77	7.11	6.8			
λ _ι (cm)	22.8	20.7	21.0	24.3	23.1			
n @ λ _{max}	2.15	2.20	1.82	2.15	2.06			
τ _{decay} (ns)	300	30/10	?	30	100			
λ _{max} (nm)	480	425/420	?	420	470			
Cut-off λ (nm)	310	350	250	280	300			
Light Output (%)	100	1.4/0.37	?	17	20			
Melting point (°C)	1050	1123	842	608	1030			
Raw Material Cost (%)	100	49	29	29	47			
IEEE Trans. Nucl. Sci. 59 (2012) 2229-2236								




Search for Scintillation in Doped PbF₂




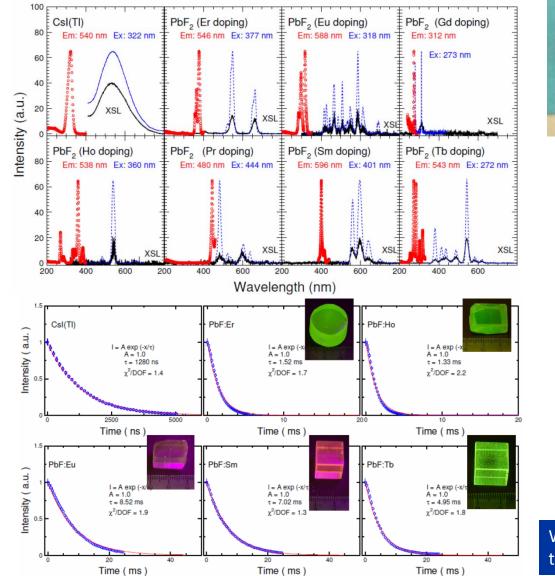

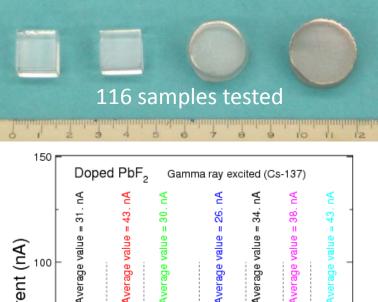

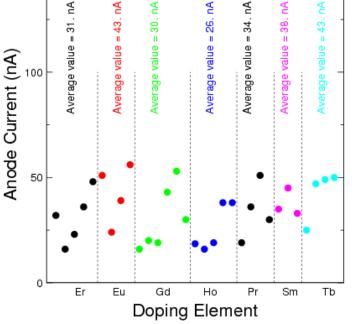

Will look performance at low temperature with the FLS920 fluorescence lifetime spectrometer

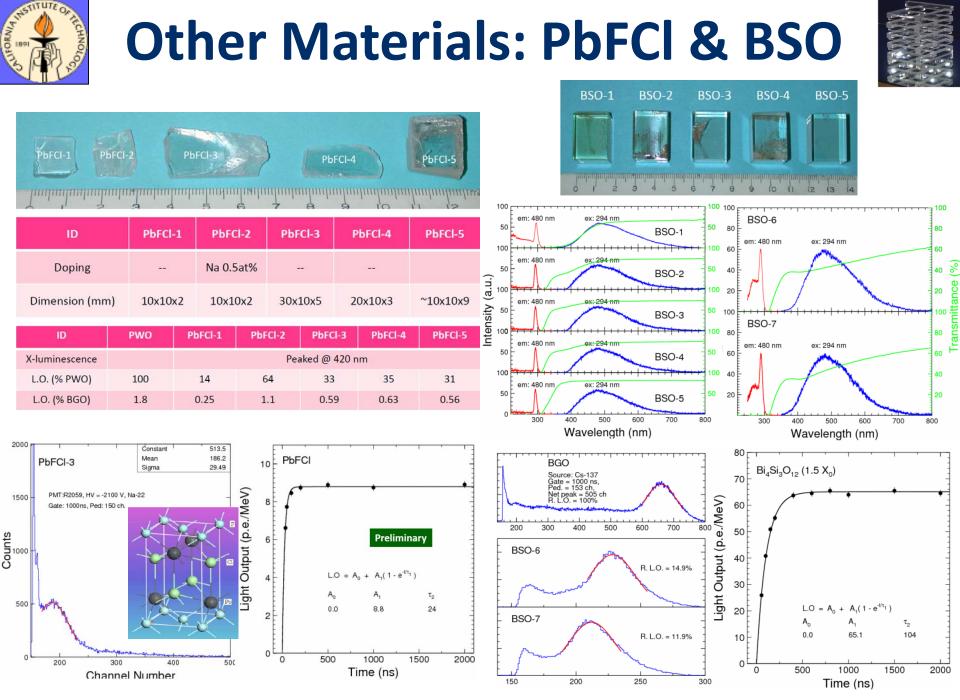

May 29,, 2015




Large Size BSO Sample





20 x 20 x 200 mm BSO crystals shows good optical and scintillation properties.

Presented in SORMA2014.

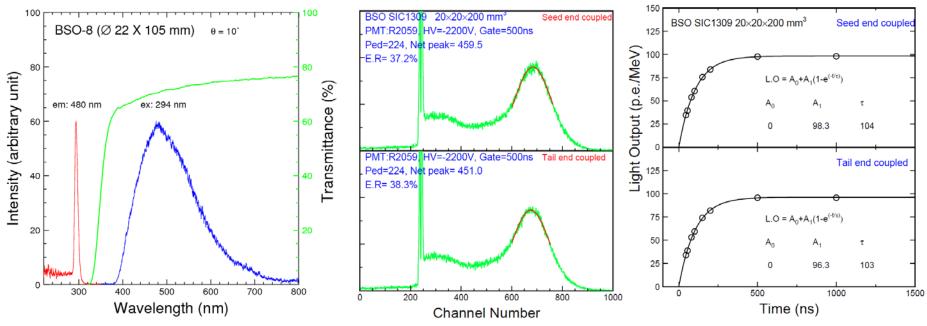



Summary



- Bright, fast and radiation hard LSO/LYSO crystals are base-lined for a total absorption ECAL for COMET at J-PARC. LYSO crystal based sampling calorimeter offers a robust stable calorimeter for the proposed HL-LHC.
- Future crystal calorimeters with more than ten times faster rate/timing capability require very fast crystals, e.g. BaF₂ with a sub-ns scintillation component.
- Cost effective crystals (PbF₂, PbFCl & BSO) may provide a foundation for a homogeneous hadron calorimeter with dual readout for both Cherenkov and scintillation light to achieve good jet mass resolution for ILC/CLIC.
- Novel scintillators in crystals, ceramics and glasses, may play important role for future HEP experiments.



LYSO Based Shashlik Cell Design



		LHCb	Plan-1	Plan-2
		Lead (Pb)	Lead (Pb)	Tungsten (W)
Absorber	Density (g/cm3)	11.4	11.4	19.3
	Radiation Length (cm)	0.56	0.56	0.35
	Moliere Radius (cm)	1.60	1.60	0.93
	dE/dX (MeV/cm)	12.74	12.74	22.1
	Thickness (mm)	2	4	2.5
	Plates number	66	28	28
		BASF-165 Polystyrene (Sc)	LYSO	LYSO
	Density (g/cm3)	1.06	7.4	7.4
	Light Yield (photons/MeV)	5200	30000	30000
Cointillaton	Radiation length (cm)	41.31	1.14	1.14
Scintillator	Moliere Radius (cm)	9.59	2.07	2.07
	dE/dX (MeV/cm)	2.05	9.55	9.55
	Plate Thickness(mm)	4	2	2
	Plates number	67	29	29
		Kurarray Y-11(250)	Kurarray Y-11(250)	Kurarray Y-11(250)
WLS Fiber	Diameter (mm)	1.2	1.2	1.2
	Number /Cell	16	4	4
	Total Depth (X0)	24.22	25.09	25.09
	Sampling Fraction (MIPs)	0.25	0.28	0.26
	Total Physical Length (cm)	40	17	12.8
	Total Sc Length (cm)	26.8	5.8	5.8
	Absorber Weight Ratio	0.84	0.75	0.76
Cell Properties	Scintillator Weight Ratio	0.16	0.25	0.24
Cell Properties	Average Density (g/cm3)	4.47	10.04	13.91
	Average Radiation Length (cm)	1.65	0.68	0.51
	Average Moliere Radius (cm)	3.6	1.7	1.2
	Transverse Dimension (cm)	4.1	1.9	1.4
	Sc-depth/Total-depth in X0	0.0268	0.2028	0.2028
	WLS Fiber Density (N/cm2)	0.97	1.06	2.07
MIPs Energy Deposition	Sc plates (MeV)	54.94	55.39	55.39
Light Yield using MIPs	Photon Electrons/GeV	3077	17897	17897
Signal of MIPs	Photon Electrons / MIP	169	991	991
Module Properties	Energy Resolution (a, %)	8.2	9.0*	9.0*

* Based on the simulation of Zhigang Wang, IHEP, Beijing.

May 29,, 2015