

Charged particle flow measured with the KRATTA detector in the ASY-EOS experiment

Sebastian Kupny, Janusz Brzychczyk, Jerzy Łukasik, Piotr Pawłowski and the ASY-EOS Collaboration

M. Smoluchowski Institute of Physics, Jagiellonian University, PI-30059 Kraków Poland
IFJ-PAN, PI-31342 Kraków Poland

sebastian.kupny@uj.edu.pl

The KRATTA detector in the ASY-EOS experiment

Sebastian Kupny^{a†}, Janusz Brzychczyk^a, Jerzy Łukasik^b, Piotr Pawłowski^b and the ASY-EOS Collaboration

^aM. Smoluchowski Institute of Physics, Jagiellonian University, PI-30059 Kraków Poland
^bIFJ-PAN, PI-31342 Kraków Poland
†sebastian.kupny@uj.edu.pl

INTRODUCTION

The KRATTA detector (Krakow Triple Telescope Array) [1] has been used to measure the energy, emission angles and isotopic composition of light charged reaction products in the ASY-EOS experiment at GSI laboratory [2].



Fig. 1 KRATTA detector during the ASY-EOS experiment (May 2013)

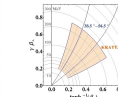


Fig. 2 Acceptance for protons and angular range used for comparison with FOPI

This versatile, low threshold, broad energy range system consists of 38 independent modules which can be arranged in an arbitrary configuration - Fig. 1 presents the used 7x5 array. A single module, covering actively about 4.5 msr of the solid angle at the optimal distance of 40 cm from the target, consists of three identical, 500 μm thick, large area photodiodes, used also for direct detection, and of two CsI(1500 ppm TI) crystals of 2.5 and 12.5 cm length, respectively. The presented results are from the analysis of the experimental data for Au+Au at 400 MeV/nucleon incident energy.

PULSE SHAPE PROCESSING

Signals from the KRATTA modules have been stored for the off-line analysis using the 100 MHz V1724 CAEN digitizers. The pulse shape analysis allowed to decompose the complex signals from the middle photodiode (Single Chip Telescope [3] segment, SCT) into the ionization and scintillation components (Fig. 3) and to obtain a satisfactory isotopic resolution with a single readout channel.

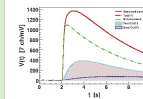


Fig. 3 Registered pulse shapes and its decomposition

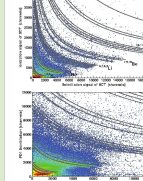


Fig. 4 Top: Identification spectra for the SCT segment. Bottom: A-E map for the thin and thick CsI(Tl) crystals

The isotopic resolution obtained using a single readout channel was found to compete very well with those obtained using the standard two channel readout (Fig. 4 top). The applied pulse shape analysis permitted also the identification of particles stopped in the first photodiode and the reduction of the identification threshold, due to the thickness of the first photodiode, by a factor of three. Thanks to the pulse shape analysis, it was also possible to obtain the ballistic deficit free amplitudes, which allowed for easy energy calibration and identification based on the reductions of the range-energy tables. For further analysis the parameters of the decomposed pulses were used in order to speed up the calculations.

REFERENCES

- [1] J. Łukasik et al., Nucl. Instr. Meth. A 709 (2013) 120128
- [2] P. Bassotto et al., Journal of Physics: Conference Series, 430, 5 (2013)
- [3] G. Pasquali et al., Nucl. Instr. Meth. A 301 (1991) 101
- [4] S. Kupny, et al. Acta Phys. Polon. B Suppl. 6, 1115-1119 (2013)
- [5] P. Bassotto et al. Phys. Lett. B, 697:471-476, (2011)
- [6] J.-Y. Ollitrault arXiv:nucl-ex/9711000v2
- [7] W. Reisdorf, FOPI Collaboration, private communication.

ANALYSIS AND RESULTS

Background recognition

In the identification maps, the isotopic lines are clearly visible, even though the background due to the secondary reactions in the long CsI(Tl) crystals is substantial (Fig. 4). In order to identify the events contributing to the background, a recognition procedure based on the self-organizing neural network has been developed [4]. This method groups together events with similar pulse shape descriptors. It allows to identify up to 90% of the background events. The efficiency of this method is presented in Fig. 5.

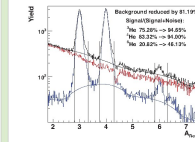


Fig. 5 Histogram isotopic distribution. Black line - raw distribution, blue - after background subtraction, red - the background.

Flow

The data measured with the KRATTA array in the ASY-EOS experiment was used to construct flow observables for light charged particles. The reaction plane orientation and the centrality of the collision were provided by the CHIMERA and ALADIN TOF Wall detectors. The reaction plane resolution corrections have been calculated according to [6]. Figure 6 shows the azimuthal distributions of protons measured with respect to the reaction plane in three rapidity regions. Target and projectile rapidity regions are dominated by the in-plane emissions, while the mid-rapidity range shows a strong squeeze-out signal.

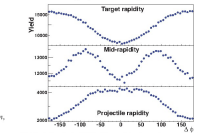


Fig. 6. Azimuthal distributions of protons measured in three regions of the rapidity: $b = (0-7.5)$ fm for Au+Au collisions at 400 MeV/nucleon.

Energy calibration

The energy calibration was performed using the characteristic punch-through points, curvatures of the $\Delta E-E$ lines, light-energy conversion formula and the ATIMA range-energy tables. Probabilities of the secondary reactions for light charged particles as a function of their range in the CsI(Tl) crystals was estimated using the Monte Carlo simulations within the GEANT4 environment. The data analysis was performed within the FairRoot framework.

Figure 7 presents the flow parameters v_n and v_n' for LCP as a function of the scaled center of mass rapidity measured for semi-central collisions with the KRATTA detector (circles) and the FOPI detector [7] (stars) in a common angular range in the CsI(Tl) crystals was estimated using the Monte Carlo simulations within the GEANT4 environment. The data analysis was performed within the FairRoot framework. The data points are corrected for the secondary reaction background. The agreement between these two detectors and experiments is very satisfactory.

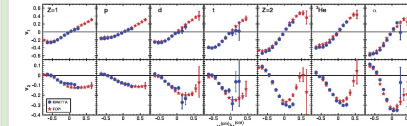


Fig. 7. Flow of LCP measured with the KRATTA (circles) and FOPI (stars) detectors. Both measurements were performed in the polar angle range $38.5^\circ < \theta < 56.5^\circ$ and for the interval of reduced impact parameter range of 0.25-0.45.

CONCLUSIONS

The KRATTA array, together with the CHIMERA and TOF-wall detectors, has proven its usefulness in flow measurements of light charged reaction products during the ASY-EOS experiment conducted at GSI. The preliminary comparison with the FOPI results measured for the same reaction in a different campaign shows a remarkable agreement and quality of the data. The results will be compared with UrQMD predictions in order to extract the information contained therein regarding the symmetry term in the nuclear equation of state.

ACKNOWLEDGMENTS

One of us (S.K.) acknowledges the support by the Foundation for Polish Science - MPD program, co-financed by the European Union within the European Regional Development Fund. Work made possible through funding by Polish Ministry of Science and Higher Education under grant No. DPN/N108/GSI/2009.



The KRATTA detecting system

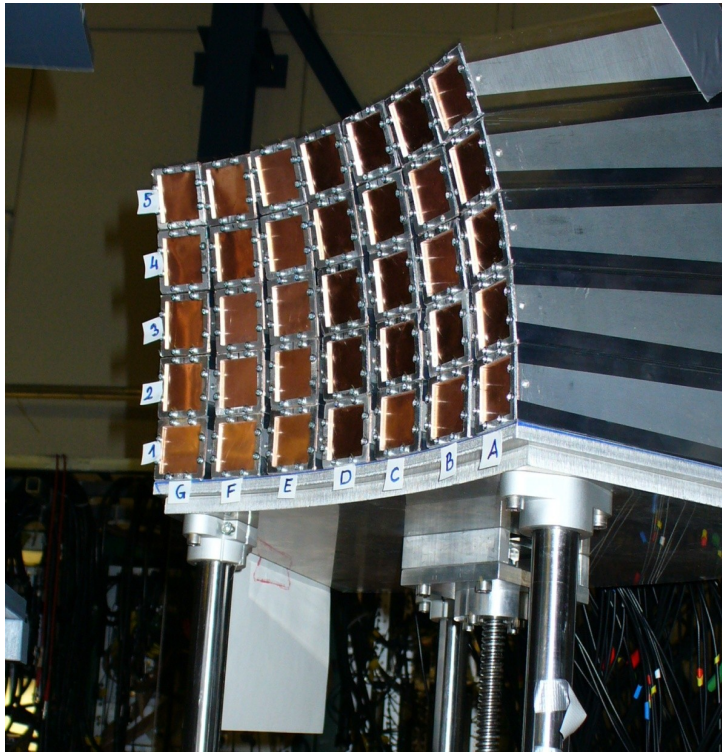


Fig.1: The KRATTA

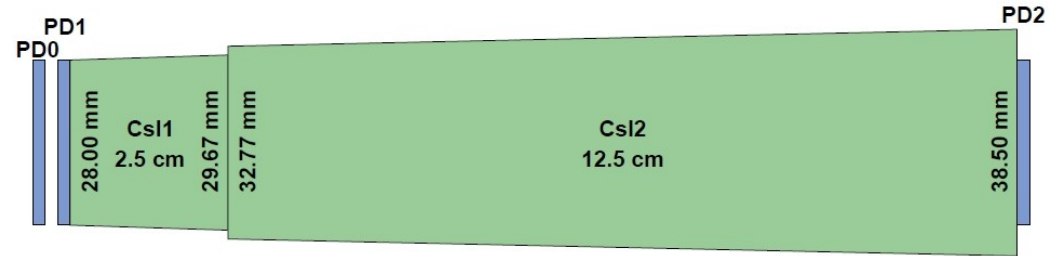


Fig.2: Active elements of the single module: 3xSilicon photodiodes and 2xCsI(Tl)

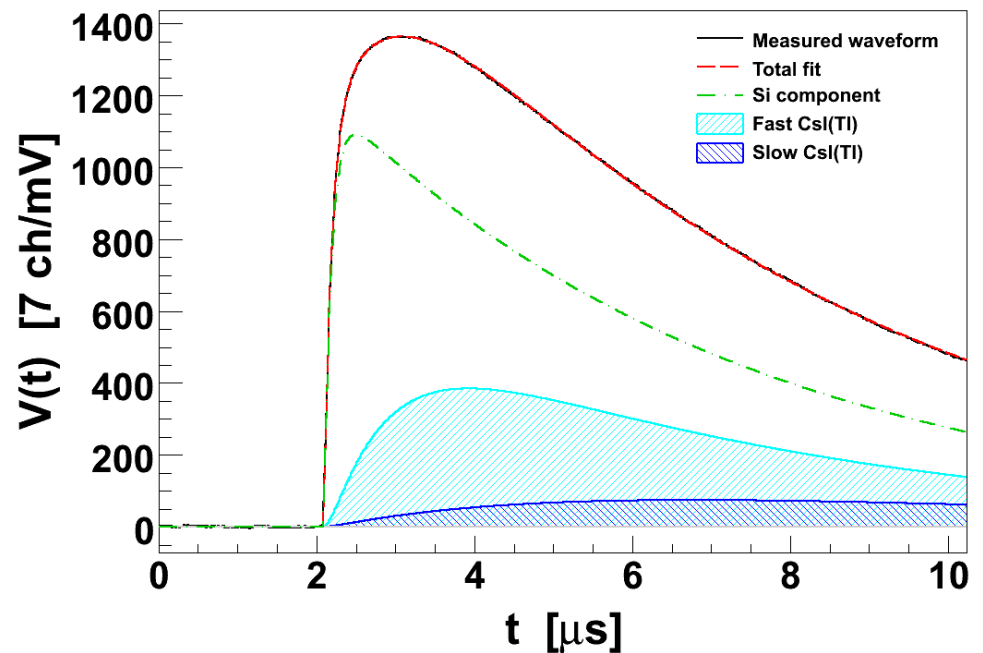


Fig.3: Pulse shape with their decomposition

The KRATTA detecting system

- Background reduction
 - With self-organizing neural networks,
 - See poster of J.Lukasik
- calibration and identification

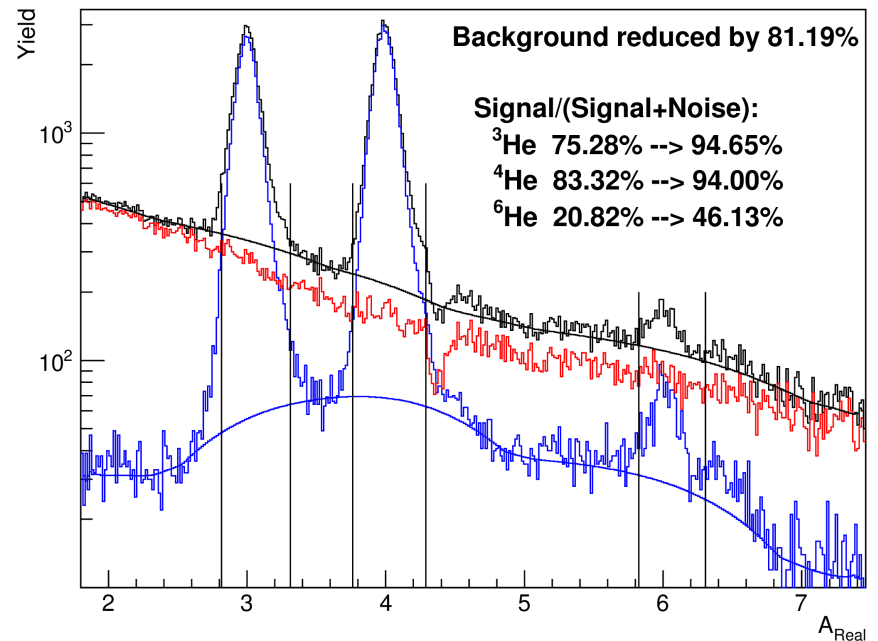


Fig.4: Helium isotopic distribution. Black line – raw distribution, blue – after background subtraction, red – the background.

Light Charged Particles flow

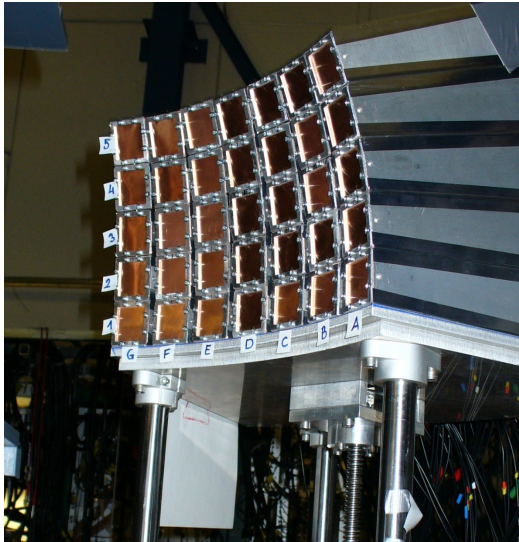


Fig.5: The KRATTA (2011)



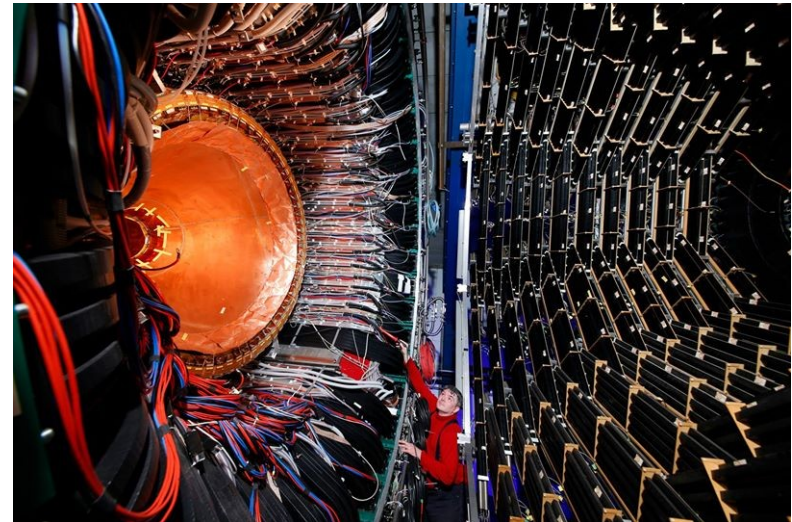


Fig.6: The FOPI detector, Foto: Thomas Ernsting/HA Hessen Agentur, 2013

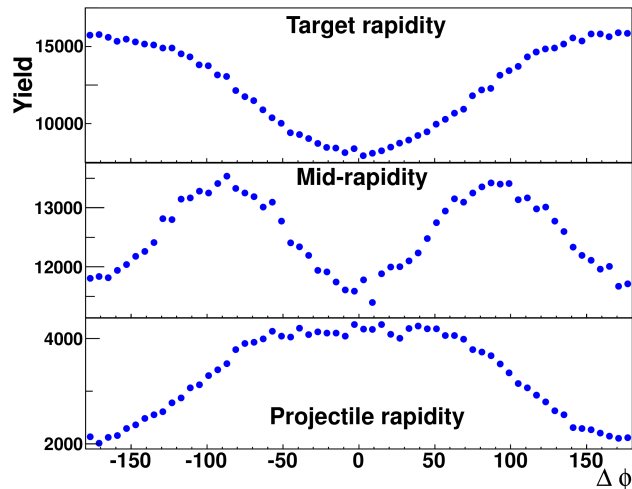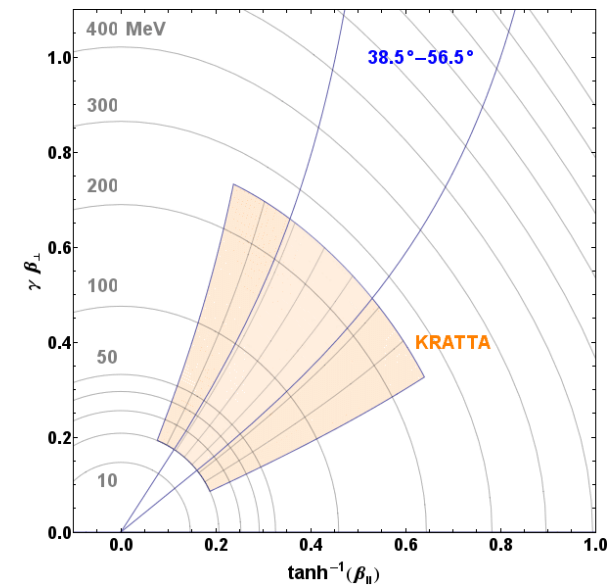

Flow of Light Charged Particles

comparison with FOPI data

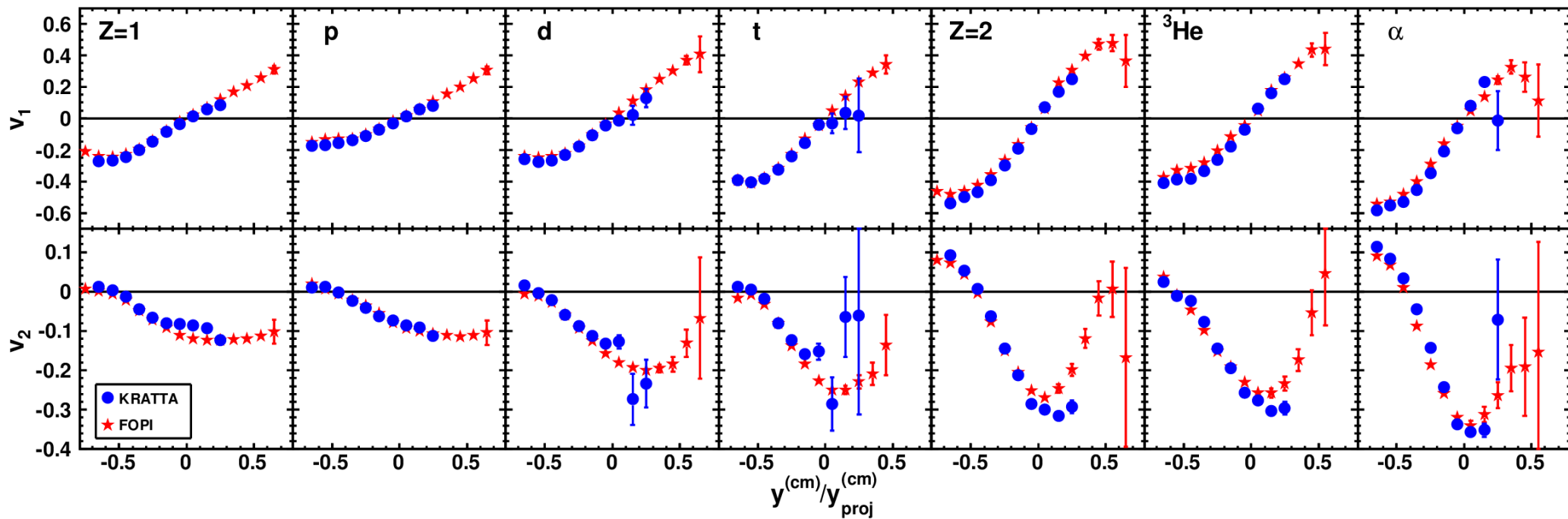


Fig.9: Flow of LCP measured with the KRATTA (circles) and FOPI (stars) detectors. Both measurements were performed in the polar angle range 38.5- 56.5 deg and for the interval of reduced impact parameter range of 0.25-0.45.